

Analysing growth curves and other user-defined plates in **opm**

Lea A.I. Vaas
Leibniz Institute DSMZ

Markus Göker
Leibniz Institute DSMZ

Abstract

This is tutorial about the analysis of growth curves and other user defined kinetics with the **opm** package in the version of December 18, 2013. It is explained how any kinds of growth or respiration measurements can be input into **opm**. We also show how Phenotype Microarray (PM) data with user-defined plate types can be analysed. Analysing such data visually and statistically requires in some cases adaptations of function arguments whose defaults are targeting PM data. All these practically relevant issues are explained in detail.

Keywords: Growth Kinetics.

1. Introduction

A detailed description of the OmniLog® Phenotype Microarray (PM) system, its measuring procedure and data characteristics are found in the vignette “**opm**: An R Package for Analysing OmniLog® Phenotype Microarray Data” (called “main tutorial” in the following). How substrate information stored within **opm** can be accessed and used for advanced visual and statistical analyses is explained in the vignette “Working with substrate information in **opm**” (called “substrate tutorial” in the following). The description of the methods below presupposes that the user is familiar with the usage of **opm** and has studied the main tutorial, the substrate tutorial as well as the entries of the **opm** manual relevant to her or his research. Especially the concepts behind, and the methods available for, the different classes of **opm** objects should be known before starting with this tutorial.

In addition to visual inspection or statistical comparative analyses of PM data, as described in the main tutorial and the substrate tutorial, users might be interested in analysing data other than PM data, or analysing PM with user-defined plate types. To work with user-defined PM plates only requires registering these plates, i.e. storing a mapping from well coordinates to substrate names, and optionally also a full, descriptive name for the plate. The analysis of data other than PM data, such as growth curves, additionally requires inputting these data and converting them to OPMX objects. Moreover, some defaults of the plotting functions are only suitable for PM data. Hence, the functions should be called slightly distinctly.

Besides these slight restrictions, which are illustrated with examples below, non-PM data can be analysed with **opm** almost like PM data.

2. Preparation

As usual, **opm** must be loaded before any analysis can begin:

```
R> if ("package:opm" %in% search())
  detach("package:opm", unload = TRUE)
R> library("opm")
```

3. Growth-curve data input

3.1. User-entered data frames

In the following we will use the growth-measurements data set from [Vaas, Marheine, Sikorski, Göker, and Schumacher \(2013\)](#) as exemplar. These data have been entered by hand and then input into R with one of the functions for reading Comma-Separated Values (CSV), yielding a data frame, which comes with **opm**:

```
R> data("potato")
R> head(potato)
```

	Genotype	Treatment	Replicate	Time	FM	DM
1	07-08-1	0.16M NaCl	1	2	597	44
2	07-08-1	0.16M NaCl	2	2	550	40
3	07-08-1	0.16M NaCl	3	2	633	48
4	07-08-1	0.16M NaCl	4	2	490	31
5	07-08-1	0.16M NaCl	5	2	617	47
6	07-08-1	0.16M NaCl	1	4	585	55

For details on this data set, enter `?potato` at the R prompt. The measurements are in “long” format and must be reshaped using the eponymous function into “wide” format. We do this separately for the Dry Mass (DM) and Fresh Mass (FM) measurements within the data set:

```
R> # order the treatment
R> potato$Treatment <- factor(potato$Treatment,
  levels = c("Control", "0.16M NaCl", "0.32M NaCl", "0.5M Sorbitol"),
  ordered = TRUE)
R> # reshape
R> potato.fm <- reshape(potato, v.names = "FM", drop = "DM", direction = "wide",
  idvar = c("Genotype", "Treatment", "Replicate"), timevar = "Time")
R> potato.dm <- reshape(potato, v.names = "DM", drop = "FM", direction = "wide",
  idvar = c("Genotype", "Treatment", "Replicate"), timevar = "Time")
```

“long” format means, that each measurement is stored in a separate record with one entry per line (see above). Thus for each data point the entries in “Genotype”, “Treatment” and “Time” have to be repeated, resulting in a data frame with dimensions of 540 rows in 6 columns.

With the `reshape` it is possible to rearrange the data set in a form, where the columns “Genotype”, “Treatment” and “Replicate” are kept and columns “Time” and either “FM”, or “DM” respectively, are merged resulting in 9 columns representing the measurement times (see below the first six rows of object `potato.fm`).

```
R> head(potato.fm)
```

	Genotype	Treatment	Replicate	FM.2	FM.4	FM.6	FM.8	FM.10	FM.12	FM.14	FM.16	FM.18
1	07-08-1	0.16M NaCl	1	597	585	882	844	1291	1847	2232	2560	2808
2	07-08-1	0.16M NaCl	2	550	614	908	1103	1240	1798	2184	2832	2501
3	07-08-1	0.16M NaCl	3	633	570	855	1200	1392	1827	2360	2522	3113
4	07-08-1	0.16M NaCl	4	490	681	1087	994	1478	1921	2315	2317	2761
5	07-08-1	0.16M NaCl	5	617	707	962	849	1446	1853	2335	2564	2426
46	07-08-1	0.32M NaCl	1	395	551	392	342	322	322	368	336	274

Thus the dimension of the data object dwindled to 60 rows in 12 columns.

Now the data are in the right arrangement for the next step, the conversion into OPMX or MOPMX objects.

The main function for converting user-defined data frames to OPMX or MOPMX objects is `opmx`, which can directly be applied to the objects created in the last step. This works because the “horizontal” input format of `opmx` corresponds to the “wide” format of `reshape`.

```
R> potato.fm <- opmx(potato.fm, position = c("Genotype", "Replicate"),
  well = "Treatment", prefix = "FM.",
  full.name = c(fm = "Growth experiment, fresh mass"))
R> potato.dm <- opmx(potato.dm, position = c("Genotype", "Replicate"),
  well = "Treatment", prefix = "DM.",
  full.name = c(dm = "Growth experiment, dry mass"))
```

The data frame contains all substrate information (in the “Treatment” column). Hence, **opm** registers the mapping from well coordinates to substrate names on the fly. The plate type must be provided, however. As it is not within the data frame, the short name of the plate type is taken from the `full.name` argument, whose main purpose is to enter the full, descriptive name of the plate type. “Genotype” and “Replicate” go to the metadata of the resulting object and together identify each plate. In the case of PM data, this is done using the position of the plate within the OmniLog® reader. Thus the relevant argument here is `position`, which must be supplied unless there is a column of that name. The `prefix` argument helps identifying the columns with measurements over time.

The registered plate type can be queried as follows:

```
R> plate_type(TRUE) # shows all existing user-defined plates
```

```
[1] "CUSTOM:DM" "CUSTOM:FM"
```

```
R> listing(wells(plate = c("CUSTOM:FM", "CUSTOM:DM")))
```

```
CUSTOM:FM:
- Growth experiment, fresh mass
- A01: 0.16M NaCl
  A02: 0.32M NaCl
  A03: 0.5M Sorbitol
  A04: Control
CUSTOM:DM:
```

- Growth experiment, dry mass
- A01: 0.16M NaCl
- A02: 0.32M NaCl
- A03: 0.5M Sorbitol
- A04: Control

Note the prefix “CUSTOM:”, which is used to distinguish user-defined plate type from those that come with **opm**. Please keep in mind, that the definition of plate types is only available in the current R session. The definitions will be lost, when the session is terminated. Saving and/or loading of a session can be managed by the functionality provided by **session** ([Warnes. 2012](#)).

The object resulting from **listing** can be output with **to_yaml** or **saveRDS** for externally storing plate types in files.

With the **potato.dm** and **potato.fm** objects the user can now follow the work flow as it is envisaged for processing of usual PM data. Please continue in Section 4 and following for plotting and statistical analysis of the estimated curve parameters.

3.2. Direct registration of plate types

An example input file comes along with **opm**, providing growth curve data which were derived from an growth challenging experiment with two *Escherichia coli* strains (Deutsche Sammlung von Mikroorganismen (DSMZ) 18039 = K12 and the type strain DSM 30083^T) on increasing Glucose concentrations. Each strain-Glucose-concentration combination was repeated twice on the plate. Thus, it will first be shown, how to prepare a plate map, register it as a new plate type and import the data with subsequently conversion of the data into an **OPMX** or **MOPMX** object. Afterwards it will be shown how to use the **split** function in order to split the objects to provide objects representing the repetitions.

3.2.1. Set up a plate map and register plate types

Provided by **register_plate**, **opm** brings several options for setting up a user defined plate layout. This function works with both customised PM plates run, as well as plates of other well design and formats used for measurement of for example growth in a conventional plate reader.

For small set ups it might be feasible to type the substrate allocation manually into a character string, as it is done in the following.

```
R> register_plate(growth = c(A01 = "negative control #1", A02 = "10mM Glucose #1",
  A03 = "20mM Glucose #1", A04 = "50mM Glucose #1", A05 = "100mM Glucose #1",
  A06 = "200mM Glucose #1",
  B01 = "negative control #2", B02 = "10mM Glucose #2",
  B03 = "20mM Glucose #2", B04 = "50mM Glucose #2", B05 = "100mM Glucose #2",
  B06 = "200mM Glucose #2",
  C01 = "negative control #3", C02 = "10mM Glucose #3",
  C03 = "20mM Glucose #3", C04 = "50mM Glucose #3", C05 = "100mM Glucose #3",
  C06 = "200mM Glucose #3",
  D01 = "negative control #4", D02 = "10mM Glucose #4",
  D03 = "20mM Glucose #4", D04 = "50mM Glucose #4", D05 = "100mM Glucose #4",
```

```

D06 = "200mM Glucose #4"),
growth = "growth")

CUSTOM:GROWTH CUSTOM:GROWTH
      TRUE      TRUE

R> #listing(wells(plate = "custom:growth"))

```

However, this procedure is error prone and not efficient when dealing with sets containing more than a few wells. Alternatively, a user-designed plate can also be registered with a plate map given as matrix. The matrix then directly represents the allocation of the used substrates on the plate.

First, the matrix for the plate map has to be stated as follows.

```

R> growth <- matrix(rep(c("negative control", "10mM Glucose", "20mM Glucose",
      "50mM Glucose", "100mM Glucose", "200mM Glucose"), each = 4),
      nrow = 4, ncol = 6)
R> rownames(growth) <- LETTERS[1:4]
R> colnames(growth) <- 1:6

```

The next step is to register the layout as a new plate type. Here it is named “growth”.

```

R> growth.reg <- register_plate(growth = growth)
R> # listing(wells(plate = "CUSTOM:GROWTH"))

```

Alternatively, the input for `register_plate` can also be a data frame.

```

R> growth.dat <- as.data.frame(growth)
R> growth.dat.reg <- register_plate(growth.dat = growth.dat)
R> # listing(wells(plate = "CUSTOM:GROWTH-DAT"))

```

3.3. Input of TECAN data

Here we will use an exemplar that comes with **opm** as input data file:

```

R> tecan.file <- opm_files("growth")
R> tecan.file <- grep("tecan", tecan.file, ignore.case = TRUE, value = TRUE)
R> tecan <- read.table(tecan.file)
R> head(tecan)

```

```

      V1      V2      V3      V4      V5      V6      V7
1 <> 1.000 2.000 3.000 4.000 5.000 6.000
2 A 0.087 0.088 0.087 0.088 0.085 0.084
3 B 0.087 0.088 0.087 0.086 0.087 0.085
4 C 0.083 0.082 0.081 0.083 0.079 0.077
5 D 0.083 0.083 0.081 0.082 0.080 0.079
6 <> 1.000 2.000 3.000 4.000 5.000 6.000

```

This format is not particularly useful within R but can be converted using the “rectangular” mode of **opmx**.

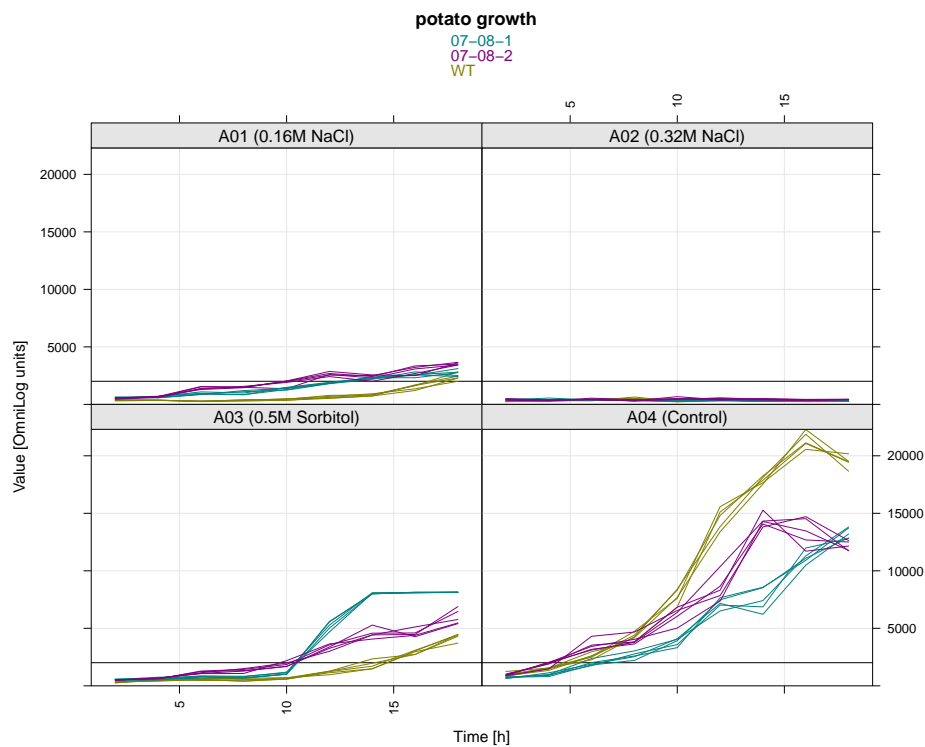


Figure 1: here some caption about the plot

```
R> # rectangular input, as it comes from the 'tecan' object
R> (y <- opmx(tecan, "rectangular", plate.type = "growth", position = 1,
  interval = 1))
```

Class	OPM
From file	
Hours measured	71
Number of wells	24
Plate type	CUSTOM:GROWTH
Position	1
Setup time	Wed Dec 18 17:14:22 2013
Metadata	0
Aggregated	FALSE
Discretized	FALSE

TODO.

4. Visualisation of growth curves

visualisation of potato data

here some description

TODO.

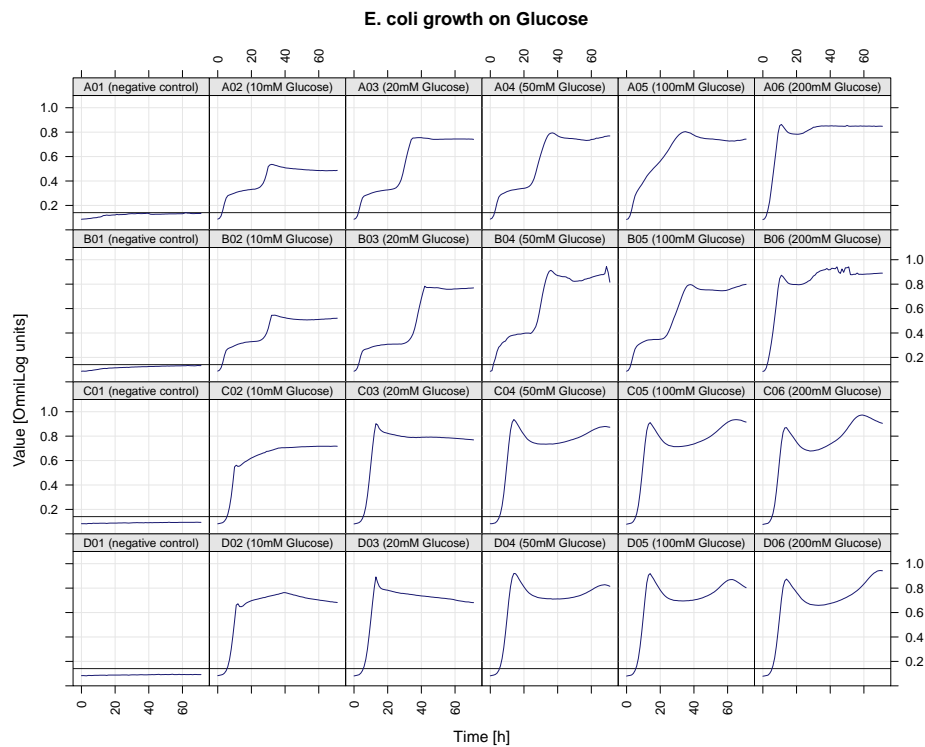


Figure 2: here some caption about the plot

```
R> print(xy_plot(potato.fm, theor.max = FALSE, include = "Genotype",
  main = list(in.parens = FALSE), neg.ctrl = FALSE, ylab = "Mass [g]"))
R> print(xy_plot(potato.dm, theor.max = FALSE, include = "Genotype",
  main = list(in.parens = FALSE), neg.ctrl = FALSE, ylab = "Mass [g]"))
```

5. Estimating parameters from growth curves

TODO.

6. Statistical analysis of growth curves

TODO.

7. Acknowledgements

We are grateful to Victoria Michael (Deutsche Sammlung von Mikroorganismen und Zellkulturen (DSMZ)) for providing growth curves measured with a TECAN instrument.

References

Vaas LAI, Marheine M, Sikorski J, Göker M, Schumacher HM (2013). “Impacts of pr-10a Overexpression at the Molecular and the Phenotypic Level.” *International Journal of Molecular Sciences*, **14**, 15141–15166. doi:10.3390/ijms140715141. URL <http://www.mdpi.com/1422-0067/14/7/15141>.

Warnes GR (2012). *session: Functions for interacting with, saving and restoring R sessions*. R package version 1.0.3, URL <http://CRAN.R-project.org/package=session>.

Affiliation:

Markus Göker

Leibniz Institute DSMZ – German Collection of Microorganisms and Cell Cultures
Braunschweig

Telephone: +49/531-2616-272

Fax: +49/531-2616-237

E-mail: markus.goeker@dsmz.de

URL: www.dsmz.de